

# DOSING PUMP FOR LIQUIDS

**SIMDOS 10**

- Simple and intuitive operation
- Easy calibration
- Adaptable to fluid characteristics
- Self priming, can run dry
- Analog and impulse control (RC)
- Splashproof IP65S
- Small footprint



## Fields of application

- Chemistry
- Synthetic organic chemistry
- Pharmacology
- Solvents for adhesives and glue
- Polymers & polymer coating
- Food technology research
- Detergents

## Hydraulic data

Flow rate	Dosing volume	Suction height	Max pressure
1 - 100 ml/min	1 - 999 ml	3 mWg	6 bar

## Chemical resistance

Code	Pump head	Diaphragm	Valves	Filter
KT	PP	PTFE covered	FFKM Kalrez®	PEEK
TT	PVDF	PTFE covered	FFKM Kalrez®	PVDF
FT	PTFE	PTFE covered	FFKM Kalrez® *	PEEK
ST	STAINLESS STEEL	PTFE covered	FFKM Kalrez®	PEEK

\*On request FFKM Chemraz®

## Functions

S Version	RC Version	RCP Version
<ul style="list-style-type: none"> <li>• Manual operation</li> </ul>	<ul style="list-style-type: none"> <li>• Manual operation</li> <li>• Analog control: 0-5 V, 0-10 V, 0-20 mA, 4-20 mA, from 1 to 100%</li> <li>• Start/Stop through logic control (TTL)</li> </ul>	<ul style="list-style-type: none"> <li>• Reset/Prime through logic control (TTL)</li> <li>• Output signal fault alarm</li> <li>• Cable for external control included</li> <li>• Same as RC</li> <li>• Plus RS 232</li> </ul>


## Technical data

- Main supply 100 - 240 V / 50 - 60 Hz
- Allowed ambient temperature +5 to +40° C
- Allowed liquid temperature +5 to +80° C
- Maximum viscosity allowed 150 cSt
- Accuracy +/-2% (nominal value)
- Repeatability +/-1%
- Power consumption 24 W
- Protection class IP65
- Dimensions: 150 x 93 x 144 mm
- Weight: 0.9 kg
- Hydraulic connection: compression fittings for 4/6 mm tube

# ACCESSORIES

## Valves Kit

	ID-Nr.
Chemraz® special valves for certain chemicals	168037

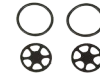

⚠ Special Valves CHEMRAZ®

**Diaphragm pressure control valve with support** when the feeding tank is above the pump or to dose into vacuum

	ID-Nr.
Kit FDV 30 KTZ	166283
Kit FDV 30 TTZ	166284
Kit FDV 30 FTZ	166285



## Filter

Typ	Material	Connector	ID-Nr.
FS 25 T	PVDF	for tubes ID 3.2 and 4 mm	165211
FS 25 X	PEEK	for tubes ID 3.2 and 4 mm	165213



**Pulsation damper** only use in continuous transfer mode from 30 to 100 ml/min

	ID-Nr.
FPD 10 KT SIMDOS	167817
FPD 10 TT SIMDOS	167818
FPD 10 FT SIMDOS	167819



## Mounting plate

	ID-Nr.
for wall mounting	160473



## Standfitting

	ID-Nr.
for retort stand	160474



## Footswitch

	ID-Nr.
for impulse start and stop	155872



SIMDOS® is a registered trademark of KNF Flodos AG. All rights reserved.  
 DuPont™ Kalrez® is a registered trademark of E. I. du Pont de Nemours and company or one of its subsidiaries. All rights reserved.  
 Chemraz® is a registered trademark of Greene, Tweed & Co. All rights reserved.

## SALES CENTERS

### Germany

KNF Neuberger GmbH  
D-79112 Freiburg  
Tel. +49 7664 5909 0  
info@knf.de  
www.knf.de

### Benelux Netherlands

KNF Verder B.V.  
NL-3451 GG Vleuten  
Tel. +31 30 677 92 40  
info@knf-verder.nl  
www.knf-verder.nl

### Benelux Belgium and Luxemburg

KNF Verder N.V.  
B-2630 Aartselaar  
Tel. +32 3 871 96 24  
info@knf.be  
www.knf.be

### China

KNF Technology (Shanghai) Co., Ltd.  
Shanghai 201203  
Tel. +86 21 5109 9695  
info@knf.com.cn  
www.knf.com.cn

### France, Morocco, Tunisia, Algeria

KNF Neuberger SAS  
F-68128 Village-Neuf  
Tel. +33 389 70 35 00  
info@knf.fr  
www.knf.fr

### UK

KNF Neuberger U.K., Ltd.  
Witney, Oxfordshire OX28 4FA  
Tel. +44 1993 77 83 73  
info@knf.co.uk  
www.knf.co.uk

### India

KNF Pumps + Systems (India) Pvt. Ltd.  
Hinjewadi, Pune 411 057  
Tel. +91 20 640 13 923  
+91 20 640 08 923  
info@knfpumps.in  
www.knfpumps.in

### Italy

KNF ITALIA S.r.l.  
I-20063 Cernusco s. Naviglio MI  
Tel. +39 02 272 03 860  
info@knf.it  
www.knf.it

### Japan

KNF Japan Co.Ltd.  
Tokyo 104-0033  
Tel. +81 3 3551 7931  
info@knf.co.jp  
www.knf.co.jp

### Korea

KNF Neuberger Ltd.  
135-502, Seoul  
Tel. +82 2 959 0255  
knf@knfkorea.com  
www.knfkorea.com

### Sweden, Finland, Denmark, Norway

KNF Neuberger AB  
SE-11743 Stockholm  
Tel. +46 8 744 51 13  
info@knf.se  
www.knf.se

### Switzerland

KNF Neuberger AG  
CH-8362 Balterswil  
Tel. +41 71 973 99 30  
knf@knf.ch  
www.knf.ch

### Taiwan

KNF Neuberger Ltd.  
Taipei City, 11490  
Tel. +886 2 2794 1011  
knftwn@knftwn.com.tw  
www.knftwn.com.tw

### USA, Canada

KNF Neuberger, Inc.  
Trenton, New Jersey 08691-1810  
Tel. +1 609 890 8600  
knfusa@knf.com  
www.knfusa.com

### Latin America

Tel. +1 609 649 1010  
gb@knf.com  
www.knf.com/es